Gene Information

Gene Symbol COL5A3 Organism Human Gene Summary This gene encodes an alpha chain for one of the low abundance fibrillar collagens. Fibrillar collagen molecules are trimers that can be composed of one or more types of alpha chains. Type V collagen is found in tissues containing type I collagen and appears to regulate the assembly of heterotypic fibers composed of both type I and type V collagen. This gene product is closely related to type XI collagen and it is possible that the collagen chains of types V and XI constitute a single collagen type with tissue-specific chain combinations. Mutations in this gene are thought to be responsible for the symptoms of a subset of patients with Ehlers-Danlos syndrome type III. Messages of several sizes can be detected in northern blots but sequence information cannot confirm the identity of the shorter messages. Gene Aliases Not Available RefSeq Accession No. NC_000019.9, NT_011295.11 UniGene ID Hs.235368 Ensembl Gene ID ENSG00000080573		
OrganismHumanGene SummaryThis gene encodes an alpha chain for one of the low abundance fibrillar collagens. Fibrillar collagen molecules are trimers that can be composed of one or more types of alpha chains. Type V collagen is found in tissues containing type I collagen and appears to regulate the assembly of heterotypic fibers composed of both type I and type V collagen. This gene product is closely related to type XI collagen and it is possible that the collagen chains of types V and XI constitute a single collagen type with tissue-specific chain combinations. Mutations in this gene are thought to be responsible for the symptoms of a subset of patients with Ehlers-Danlos syndrome type III. Messages of several sizes can be detected in northern blots but sequence information cannot confirm the identity of the shorter messages.Gene AliasesNot AvailableRefSeq Accession No.NC_000019.9, NT_011295.11UniGene IDHs.235368Ensembl Gene IDENSG00000080573	Gene Name	collagen, type V, alpha 3
Gene Summary This gene encodes an alpha chain for one of the low abundance fibrillar collagens. Fibrillar collagen molecules are trimers that can be composed of one or more types of alpha chains. Type V collagen is found in tissues containing type I collagen and appears to regulate the assembly of heterotypic fibers composed of both type I and type V collagen. This gene product is closely related to type XI collagen and it is possible that the collagen chains of types V and XI constitute a single collagen type with tissue-specific chain combinations. Mutations in this gene are thought to be responsible for the symptoms of a subset of patients with Ehlers-Danlos syndrome type III. Messages of several sizes can be detected in northern blots but sequence information cannot confirm the identity of the shorter messages. Gene Aliases Not Available RefSeq Accession No. NC_000019.9, NT_011295.11 UniGene ID Hs.235368 Ensembl Gene ID ENSG00000080573	Gene Symbol	COL5A3
Fibrillar collagen molecules are trimers that can be composed of one or more types of alpha chains. Type V collagen is found in tissues containing type I collagen and appears to regulate the assembly of heterotypic fibers composed of both type I and type V collagen. This gene product is closely related to type XI collagen and it is possible that the collagen chains of types V and XI constitute a single collagen type with tissue-specific chain combinations. Mutations in this gene are thought to be responsible for the symptoms of a subset of patients with Ehlers-Danlos syndrome type III. Messages of several sizes can be detected in northern blots but sequence information cannot confirm the identity of the shorter messages.Gene AliasesNot AvailableRefSeq Accession No.NC_000019.9, NT_011295.11UniGene IDHs.235368Ensembl Gene IDENSG0000080573	Organism	Human
RefSeq Accession No. NC_000019.9, NT_011295.11 UniGene ID Hs.235368 Ensembl Gene ID ENSG0000080573	Gene Summary	Fibrillar collagen molecules are trimers that can be composed of one or more types of alpha chains. Type V collagen is found in tissues containing type I collagen and appears to regulate the assembly of heterotypic fibers composed of both type I and type V collagen. This gene product is closely related to type XI collagen and it is possible that the collagen chains of types V and XI constitute a single collagen type with tissue-specific chain combinations. Mutations in this gene are thought to be responsible for the symptoms of a subset of patients with Ehlers-Danlos syndrome type III. Messages of several sizes can be detected in northern blots but sequence
UniGene ID Hs.235368 Ensembl Gene ID ENSG0000080573	Gene Aliases	Not Available
Ensembl Gene ID ENSG0000080573	RefSeq Accession No.	NC_000019.9, NT_011295.11
	UniGene ID	Hs.235368
Entrez Gene ID 50509	Ensembl Gene ID	ENSG0000080573
	Entrez Gene ID	50509

Assay Information

, , , , , , , , , , , , , , , , , , ,		
Unique Assay ID	qHsaCEP0058149	
Assay Type	Probe - Validation information is for the primer pair using SYBR [®] Green detection	
Detected Coding Transcript(s)	ENST0000264828	
Amplicon Context Sequence	CCCTCTGTAGTCCGACCTCAGTCCTTTCACCGTGGAAGTAGAAAGGATCAAAGG GGGTGGGGGTAGGGGTCCCAGCCTGGGATCAGGAGACATCCAGTGCTTAACTC TGGGGAGGGGTAGCACCAAAAGGTGTGAA	
Amplicon Length (bp)	106	
Chromosome Location	19:10070578-10070713	
Assay Design	Exonic	
Purification	Desalted	

Validation Results

Efficiency (%)	104
R ²	0.9976



cDNA Cq	24.53
cDNA Tm (Celsius)	84.5
gDNA Cq	24.32
Specificity (%)	100


Information to assist with data interpretation is provided at the end of this report.



COL5A3, Human


Amplification Plot

Amplification of cDNA generated from 25 ng of universal reference RNA

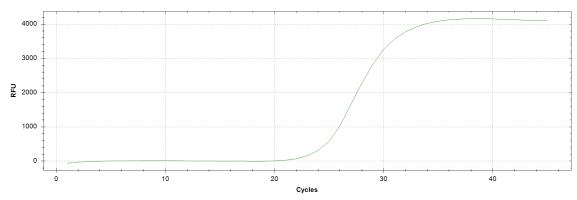


Melt Peak

Melt curve analysis of above amplification

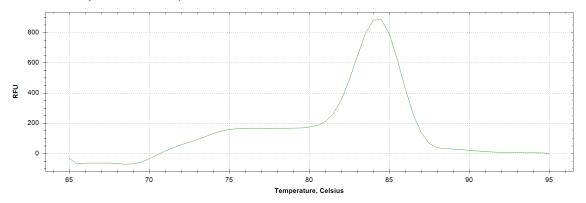

Standard Curve



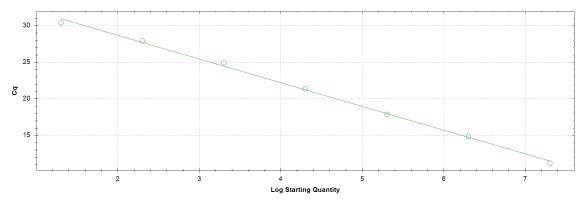



Real-Time PCR Instrument	CFX384 Real-Time PCR Detection System
Reverse Transcription Reagent	iScript™ Advanced cDNA Synthesis Kit for RT-qPCR
Real-Time PCR Supermix	SsoAdvanced [™] SYBR® Green Supermix
Experimental Sample	qPCR Human Reference Total RNA

Products used to generate validation data

Data Interpretation

-	
Unique Assay ID	This is a unique identifier that can be used to identify the assay in the literature and online.
Detected Coding Transcript(s)	This is a list of the Ensembl transcript ID(s) that this assay will detect. Details for each transcript can be found on the Ensembl website at www.ensembl.org.
Amplicon Context Sequence	This is the amplicon sequence with additional base pairs added to the beginning and/or end of the sequence. This is in accordance with the minimum information for the publication of real-time quantitative PCR experiments (MIQE) guidelines. For details, please refer to the following publication, "Primer Sequence Disclosure: A Clarification of the MIQE Guidelines" (Bustin et al 2011).
Chromosome Location	This is the chromosomal location of the amplicon context sequence within the genome.
Assay Design	Exonic: Primers sit within the same exon in the mRNA transcript and can potentially co-amplify genomic DNA. If performing gene expression analysis, it is suggested that the samples be treated with a DNase to eliminate potential unwanted signal from contaminating genomic DNA.
	Exon-exon junction: One primer sits on an exon-exon junction in mRNA. When performing gene expression analysis, this design approach will prevent unwanted signal from contaminating genomic DNA.
	Intron-spanning: Primers sit within different exons while spanning a large intron in the mRNA (intron is greater than 750bp). When performing gene expression analysis, this design approach should limit potential unwanted signal from contaminating genomic DNA.
	Small intron-spanning: Primers sit within different exons with a short intron in between (intron is smaller than 750bp). Small introns may not prevent unwanted signal from contaminating genomic DNA.
Efficiency	Assay efficiency was determined using a seven-point standard curve from 20 copies to 20 million copies. While an efficiency of 100% represents a perfect doubling of template at every cycle and is ideal, typical ranges of good assay efficiency are between 90-110%. For difficult targets, assay efficiency outside of this range are accepted and reported accordingly.
R ²	The R ² represents the linearity of the standard curve and how well the standard curve data points fit the linear regression line. Acceptable values are >0.98.



cDNA Cq	Cq value obtained from 25ng of cDNA transcribed from universal RNA when performing wet-lab validation of the assay.
	Note: Not all genes will be expressed at a detectable level in the universal RNA sample.
cDNA Tm	Melting temperature of the amplicon when running a melt curve analysis.
gDNA Cq	Cq value obtained when running the assay with 2.5ng of genomic DNA. This is more than a moderate level of genomic DNA contamination. Intron-spanning and exon-exon junction assay designs can minimize or eliminate genomic DNA detection.
	Note: Genomic DNA contamination is often present at variable levels. If concerned about genomic DNA contamination, the genomic DNA contamination control assay is recommended to run with your sample to determine if genomic DNA levels are sufficient to negatively impact results.
Specificity	This value is the percent of specific amplicon reads as measured by next generation sequencing (NGS). While 100% specificity is desirable, small decreases in specificity (<1%) can be due to NGS read errors. More significant reductions are likely due to co-amplification of homologous regions.
	Note: Since gene expression can be cell type and condition specific, the exact level and impact of co-amplification in a given sample is impossible to predict. If co-amplification is detected, it should be taken into consideration and reported when analyzing gene expression results.

